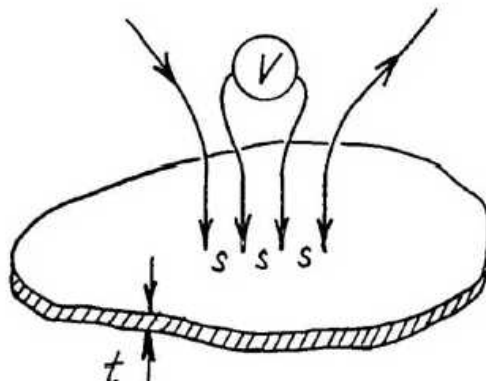


D.2) Thin Sample.

When the sample is thin, it is more convenient to write the resistivity in the form given by Smits (e) :

$$\rho = G \frac{V}{I}, \text{ where}$$

$$G = \frac{\pi}{\ln 2} \cdot t \cdot T_2\left(\frac{t}{s}\right) = 4.5324 \cdot t \cdot T_2\left(\frac{t}{s}\right); \quad (11)$$

$\frac{\pi}{\ln 2} \cdot t = 4.5324 \cdot t$ is the geometric factor for an infinitely large, thin slice (e). Thin slice means $\frac{t}{s} \ll 1$ (in practice $\frac{t}{s} \leq 0.5$).

T_2 is an additional correction factor to apply when t is not much less than s .

$$T_2 \rightarrow 1 \text{ as } \frac{t}{s} \rightarrow 0.$$

$T_2\left(\frac{t}{s}\right)$ is tabulated and plotted on the following page.